

Forma Model 900 Series

-86C ULT Upright Freezer

Operating and Maintenance Manual 7090902 Rev. 18

Visit us online to register your warranty www.thermoscientific.com/labwarranty



Important installer and user information:

A redundant temperature sensing device has been included in this ULT freezer. This device is a type "T" thermocouple. For convenient access, the thermocouple (Figure 1-3) terminates in an interconnect jack (Figure 1-5) behind the base front cover. (May be located differently in chests. See Section 1.) It is strongly recommended that this thermocouple be attached to a redundant 24 hour 7 day monitoring system with alarm capabilities. Connecting the sensor to a monitoring and alarm system separate from the freezer provides the utmost in product safety, should the integral system fail. \blacktriangle

Models Covered

Catalog No.	Capacity (cu ft)	Voltage
902	13	230/50
903	13	120/60
909	13	208-230/60 Hz
904	17	120/60
905	17	230/50
910	17	208-230/60 Hz
906	23	230/50
956	23	120/60
957	23	208-230/60 Hz
907	28	230/50
908	28	208-230/60 Hz
Double Door	(DD) Units	
991	13	230/50
992	13	120/60
986	13	208-230/60 Hz
993	17	120/60
994	17	230/50
988	17	208-230/60 Hz
990	23	120/60
995	23	230/50
989	23	208-230/60 Hz

Packing	List
---------	------

Part No.	Description	Quantity
122005	Кеу	2 (4 for DD)
380520	Neoprene Cap	2
510016	¼-20 x 5-½" Bolt	2
195763	Retaining Clip	1
370563	Remote Alarm Connector	1

MANUAL NUMBER 7090902

18	40634	10/31/17	Updated exploded drawings	bpg
17	41728	10/02/17	Added Door operation info (1-6)	bpg
16	31514	8/25/17	Updated 195964 High/Low temp stage drawings	bpg
15	41370	8/14/17	Updated BUS board information	bpg
14	41159/41348	5/12/17	Corrected oerational temp range and D-volt rating	bpg



Important Read this instruction manual. Failure to read, understand and follow the instructions in this manual may result in damage to the unit, injury to operating personnel, and poor equipment performance.

Caution All internal adjustments and maintenance must be performed by qualified service personnel. ▲

Material in this manual is for information purposes only. The contents and the product it describes are subject to change without notice. Thermo Fisher Scientific makes no representations or warranties with respect to this manual. In no event shall Thermo be held liable for any damages, direct or incidental, arising out of or related to the use of this manual.

The -86C Freezers (see page i for specific list of models) described in this manual are high performance units which can be used for research and in situations that directly support medical applications. When these products are used to support a medical application, it is an accessory to a medical device and is therefore considered as a medical device in its own right by the regulatory body (e.g. FDA).

This product is intended for use:


- As cold storage in research use.
- As a medical device for diagnostic use (storage of samples not intended to be re-introduced to human body).

Registration: This medical application is considered a Class I medical device by the FDA. This product is classified as product code – JRM, regulation number 862.2050 and is considered a Class 1 device, 510(K) exempt.

©2014 Thermo Fisher Scientific. All rights reserved.



Important operating and/or maintenance instructions. Read the accompanying text carefully.

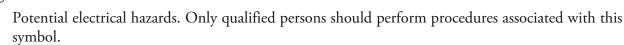



Equipment being maintained or serviced must be turned off and locked off to prevent possible injury.



Extreme temperature hazards, hot or cold. Use special handling equipment or wear special, protective clothing.



WEEE Compliance: Thermo Fisher Scientific has contracted with companies for recycling/disposal in each EU Member State. For further information, send an email to weee.recycle@thermofisher.com.

- ✓ Always use the proper protective equipment (clothing, gloves, goggles, etc.)
- ✓ Always dissipate extreme cold or heat and wear protective clothing.
- ✔ Always follow good hygiene practices.
- ✓ Each individual is responsible for his or her own safety.

Do You Need Information or Assistance on **Thermo Scientific Products?** If you do, please contact us 8:00 a.m. to 6:00 p.m. (Eastern Time) at: 1-740-373-4763 Direct 1-800-438-4851 Toll Free, U.S. and Canada 1-877-213-8051 FAX http://www.thermofisher.com Internet Worldwide Web Home Page service.led.marietta@thermofisher.com Tech Support Email Address www.unitylabservices.com Certified Service Web Page Our Sales Support staff can provide information on pricing and give you quotations. We can take your order and provide delivery information on major equipment items or make arrangements to have your local sales representative contact you. Our products are listed on the Internet and we can be contacted through our Internet home page. Our Service Support staff can supply technical information about proper setup, operation or troubleshooting of your equipment. We can fill your needs for spare or replacement parts or provide you with on-site service. We can also provide you with a quotation on our Extended Warranty for your Thermo Scientific products. Whatever Thermo Scientific products you need or use, we will be happy to discuss your applications. If you are experiencing technical problems, working together, we will help you locate the problem and, chances are, correct it yourself...over the telephone without a service call. When more extensive service is necessary, we will assist you with direct factory trained technicians or a qualified service organization for on-the-spot repair. If your service need is covered by the warranty, we will arrange for the unit to be repaired at our expense and to your satisfaction. Regardless of your needs, our professional telephone technicians are available to assist you Monday through Friday from 8:00 a.m. to 6:00 p.m. Eastern Time. Please contact us by telephone or fax. If you wish to write, our mailing address is:

Thermo Fisher Scientific (Asheville) LLC 401 Millcreek Road, Box 649 Marietta, OH 45750

International customers, please contact your local Thermo Scientific distributor.

Table of Contents

Section 1	Installation and Start-up	1-1
	Control Panel Keys, Displays, Indicators	
	Operation of the Keypad	
	Install Freezer	
	Choose Location	
	Door Operation	
	Install Wall Bumpers	
	Install Shelves	
	RS-232 Communications	1-6
	Remote Alarm Contacts	1-7
	Attach Power Cord	1-8
	Connect Unit to Electrical Power	1-8
	Run Mode	1-8
	Freezer Start-Up	1-9
	Set Operating Temperature	1-9
	Set High Temperature Alarm	1-10
	Set Low Temperature Alarm	1-10
Section 2	Calibrate	2-1
	Calibrate Control Probe	
Section 3	Alarms	3-1
	High Stage System Failure	
	Probe Failure Alarm	3-2
	Alarm Error Charts	3-3
Section 4	Maintenance	4-1
	Clean Cabinet Exterior	4-1
	Clean Air Filter	
	Clean Condenser	
	Clean Water-cooled Condenser	
	Prepare Unit for Storage	
	Long Term Storage with Water-cooled Condenser	
	Defrost Chamber	

	Clean Door Gasket	
	Vacuum Relief Port	
	Check Battery(s)	
	Replace Battery(s)	
	PREVENTIVE MAINTENANCE	
Section 5	Factory Installed Options	5-1
	Back Up System (BUS)	
	Install Vent Stack, Solenoid and Injection Assembly	
	Install Temperature Probe	
	Connect Probe/Solenoid Harness	
	BUS Control Panel	5-5
	Configure Optional BUS (Back-Up System)	
	Set Optional BUS Set Point	5-8
	Test BUS Operation	5-8
	Clean Vent Stack	
	Disconnect Fitting Asm. and Transfer Hose	5-9
	Chart Recorder	5-9
	Installing the Chart Paper	
	Change Recorder Temperature Range	
	Recorder Calibration	
	Water-cooled Condenser	
	Water Connections	
	Five Inner Door Option	
Section 6	Specifications	6-1
o o o contra o	Single Door Units	
	Double Door Units	
Section 7	Spare Parts	7 - 1
0 (* 0		
Section 8	Refrigeration Schematics	8 - 1
Section 9	Electrical Schematics	9 - 1
Section 10	Warranty	10 - 1
Appendix	Handling Liquid Nitrogen	Δ-1
- The survey	Handling Liquid CO ₂	
	First Aid	U - I

Section 1 Installation and Start-up

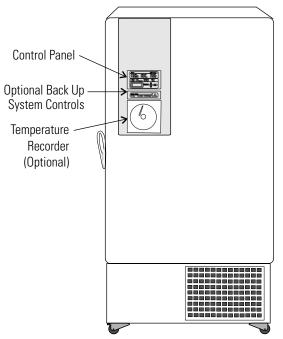


Figure 1-1. Freezer Front View Components

- Control panel keypad, displays and indicators
- BUS (Optional Back Up System) control panel
- Optional temperature recorder or datalogger

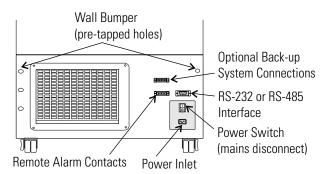


Figure 1-2. Freezer Rear View Components

- Remote alarm contacts
- Power inlet for power cord connection
- Optional BUS connections for probe and solenoid
- RS-232 (default) or RS-485 interface
- Power switch (mains disconnect)

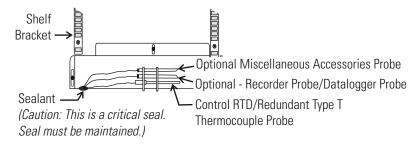


Figure 1-3. Chamber Probe(s)

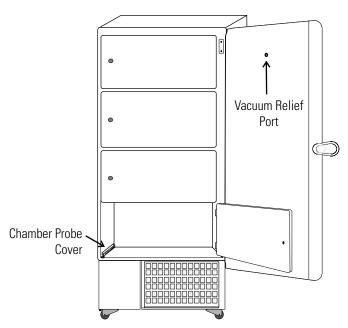


Figure 1-4. Vacuum Relief and Probe Cover Location

- Vacuum relief port pressure equalization port
- Probe cover houses control, optional recorder, datalogger, or optional miscellaneous accessory probes

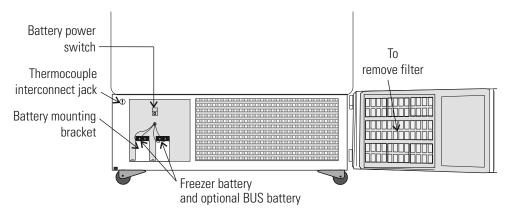


Figure 1-5. Battery(s) location and switch

- Battery mounting bracket(s)
- Battery power switch (freezer and BUS)
- Freezer battery
- Optional BUS battery
- Freezer filter location

Control Panel Keys, Displays, Indicators

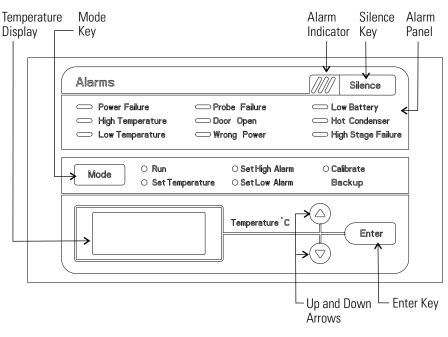


Figure 1-6. Control Panel

- Temperature Display Displays temperature in degrees Celsius.
- Mode Select Switch Used to select Run, Set Temperature, Set High Alarm, Set Low Alarm, Calibrate, Backup.
- Alarm Indicator Light pulses on/off during an alarm condition of the cabinet.
- Silence Silences the audible alarm. See Section 4 for alarm ringback times.
- Alarm Panel indicates the current alarm condition.
- Up and Down Arrows Increases or decreases values, toggles between choices.
- Enter Stores the value into memory

Operation of the Keypad

The Model 900 Series freezer has five basic modes which allow freezer setup and operation. Press the Mode key to scroll through the mode selections.



Up Arrow: Increases or toggles the parameter value.

Enter: Must press Enter key to save to memory all changed values.

Down Arrow: Decreases or toggles the parameter value.



Silence Key: Press to silence the audible alarm. See Section 4 for alarm ringback times.

Install Freezer

Note If tipped more than 45°, allow the unit to sit upright for 24 hours before start up. \blacktriangle

To remove the freezer from the pallet, use a 1/2" wrench to remove all the bolts securing the shipping bracket to the pallet.

Remove the shipping bracket. Remove the ramp boards from the pallet and place the slotted end over the ramp brackets on the pallet. The support blocks on the ramps will be facing down. Before moving the freezer, make sure the casters are unlocked and moving freely. Align the caster with the ramp boards. Use adequate personnel to roll the freezer off the pallet.

Choose Location The freezer can be easily pushed to the desired approved location, described below. If necessary, the doors and lower front panel may be opened to move the unit through tight openings. When the freezer is in position, set the front caster brakes.

Note The freezer must not be moved with the product load inside. \blacktriangle

Note For proper ventilation and airflow, a minimum clearance of 5" at the rear and top and a clearance of 8" on the side of the freezer is required. Allow adequate space in the front of the freezer for door opening. \blacktriangle

Locate the freezer on a firm, level surface in an area with an ambient temperature between 18°C and 32°C. Provide ample room to reach the mains disconnect switch (power switch) located on the rear of the freezer.

Door Operation			e equipped with an adva ltra-low temperature fre	-
	Features include:			
	• One-ha	nd operation		
	• A front-accessible lock			
	• Hasps for a standard padlock to provide additional security. Length of the shackle must be between ³ / ₄ inch (1.9 cm) and 1 ¹ / ₂ inch (3.8 cm).			
	• Durable	construction	for reliable operation ar	nd safe product storage.
Install Wall Bumpers	The parts bag, located inside the cabinet, contains the following parts.			s the following parts.
	Quantity	Stock #	Description	Purpose
	2	510016	1/4-20 x 5-1/2" Bolt	Wall Bumper
	2 Table 1-1. P	380520 Parts bag	Neoprene Cap	Cap Protector
Install Shelves	Install the s shelf level. I	nstall the shel		-
		Ũ	e optional 5 inner door o g the inner door kit. ▲	option, refer to the
RS-232 Communications	data commu	900 Series fre inications inte ult setting is F	erface. The	
	interface is s nine pin, su located on t	identification shown in Figu b "D" style co he back of the for the location	re 1-7. One Figure 1-	0 0 0_{6}
	standard DI	B9 serial exten	perature information ev sion cable can be used t devices may require a nu	o connect the freezer to a

RS-232 Communications (continued)

Data format:

Baud
Data bits 8 (7 bit ASCII with leading zero)
Start bits
Stop bits
Parity none

The data transfer sequence is transmitted in the following format. X refers to numerical temperature data.

(NUL) (-) XXX (SP) C (SP) (Error Message) (SP) (LF) (CR) (EOT) (SP)

In the event of a CNTRLFAIL, Er07, or the control probe is out of range error, the numerical temperature data (XXX) in the transmission would be replaced by T_ERR.

If no alarm condition exists, spaces will be sent. A total of 20 characters will be sent.

SP - Space	LF - Line feed
CR - Carriage return	EOT - End of text (4)
NUL - Null character (00)	

If an alarm condition does exist, "Error Message" in the protocol will be replaced by the following:

See Figure 1-2 for the location of the remote alarm contacts. The remote

UNDERTEMP (temperature above the low alarm setpoint) OVERTEMP (temperature below the high alarm setpoint) PWRFAIL (AC power failure) CNTRLFAIL (Control probe failure) Er07 (micro failure) HSHX FAIL (Heat exchanger failure) HOT COND (Hot condenser)

Remote Alarm Contacts

IMPORTANT USER INFORMATION

CAUTION! Stored product should be protected by a redundant 24 hour/day monitoring system with alarm capability. An interconnect jack and thermocouple are installed for centralized monitoring, should on-board system fail.

alarm connector is located in the parts bag provided with the manual. It must be installed if connecting the freezer to an alarm system. After installing the wiring from the alarm system to the connector, install the connector to the freezer microboard and secure with the two screws provided. The remote alarm provides a NO (normally open) output, a NC normally closed) output and COM (common). The contacts will trip on a

REMOTE CONTACTS/ANALOG OUTPUT
PIN# 1 Not Connected
PIN# 2 Not Connected
PIN# 3 Not Connected
PIN# 4 Not Connected
PIN# 5 Normally Closed
PIN# 6 Common.
PIN# 7 Normally Open
CONTACT RATING: 1A @ 30V
CONTACTS IN ALARM STATE

Figure 1-8. Remote Alarm Contacts

Remote Alarm Contacts (cont) power outage and high temperature or low temperature alarms. They will also trip on high stage, control probe or microboard failures.

Figure 1-8 shows the remote contacts in the alarm state.

Attach Power Cord

Insert the power cord into the power inlet module. Place the retaining bracket (P/N 195763) over the connector. Tighten retaining screws to secure.

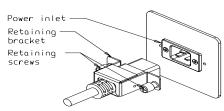


Figure 1-9. Power Cord Connection

Connect Unit to Electrical Power



Caution See the serial tag on the side of the unit for electrical specifications or refer to the electrical schematics in this manual. \blacktriangle

Caution If the factory installed option water-cooled condenser is present, do not turn the freezer on without water connected and flowing. Damage to the refrigeration system could occur within 5 minutes if water is not connected and flowing on unit start-up. Refer to Section 5. \blacktriangle

The freezer should be operated on a dedicated grounded service. Check the voltage rating on the serial tag of the unit and compare it with the outlet voltage. Then, with the power switch turned off, plug the line cord into the wall outlet.

First, turn on the freezer power switch. Then open the lower front door by grasping the bottom left corner. Locate the battery switch (Figure 1-5) and turn it to Standby mode (\circlearrowright). During initial freezer start-up, the system battery may require charging and the Low Battery indicator may illuminate.



Caution Ensure the battery switch is turned to Standby mode (\bigcirc). The rechargeable batteries require 36 hours to charge at initial start-up. A "Low Battery" alarm may occur until the batteries are fully charged. Should a power failure occur during the initial start-up period, the electronics will have limited operation.

Run Mode

Run mode is the default mode for the freezer. This mode displays the cabinet temperature on the temperature display under normal operating conditions. In addition, the Run mode allows display of the high stage heat exchange temperature.

Run Mode (cont)

Freezer Start-Up

This information scrolls by pressing the up or down arrow keys. The display returns to the operating temperature in 10 seconds if no keys are pressed.

With the freezer properly installed and connected to power, system set points can be entered. The following set points can be entered in Settings mode: Control temperature, high temperature alarm set point, low temperature alarm set point, and (optional) BUS set point. Default settings are shown in the table below.

Table 1-2. Default Settings

Default Settings	Temperature
Control Set Point	-80°C
High Temperature Alarm	-70°C
Low Temperature Alarm	-90°C
Optional BUS Set Point	-60°C

Note If the set point is changed and the low temp and high temp alarms are set 10° from the set point, the alarm set points will be adjusted automatically to maintain a distance of at least 10° from set point.

Caution If the factory installed option water-cooled condenser is present, do not turn the freezer on without water connected and flowing. Damage to the refrigeration system could occur within 5 minutes if water is not connected and flowing on unit start-up. Refer to Section 5. ▲

Set Operating Temperature Model 900 and 5900 Series freezers have an operating temperature range of -50°C to -86°C, depending on ambient temperature. The freezer is shipped from the factory with a temperature set point of -80°C. To change the operating temperature set point:

- 1. Press the Mode key until the Set Temperature indicator lights.
- 2. Press the up/down arrow key until the desired temperature set point is displayed.
- 3. Press Enter to save the set point.
- 4. Press the Mode key until the Run indicator lights for Run mode

If no keys are pressed, the freezer will automatically return to RUN mode after 5 minutes.

Note If the set point is changed and the low temperature and high temperature alarms are set 10° from the set point, the alarm set points will be adjusted automatically to maintain a distance of at least 10° from set point.

Set High Temperature Alarm	The high temperature alarm will activate an audible/visual warning when the freezer chamber temperature has reached or exceeded the high temperature alarm set point.
	To set the high temperature alarm set point:
	1. Press the Mode key until the Set High Alarm indicator lights.
	2. Press the up or down arrow key until the desired high temperature alarm set point is displayed.
	3. Press Enter to save the setting.
	4. Press the Mode key until the Run indicator lights for Run mode
	If no control keys are pressed, the freezer will automatically return to RUN mode after 5 minutes.
	Note The high alarm set point must be set at least 5°C from the control set point. \blacktriangle
	Note At initial start-up, the high temperature alarm is disabled until the cabinet reaches set point or 12 hours elapse. ▲
Set Low Temperature Alarm	The low temperature alarm will activate an audible/visual warning when the freezer chamber temperature has reached or decrease below the low temperature alarm set point.
	To set the low temperature alarm set point:
	1. Press the Mode key until the Set Low Alarm indicator lights.
	 Press the up or down arrow key until the desired low temperature alarm set point is displayed.
	alarm set point is displayed.
	alarm set point is displayed.3. Press Enter to save the setting.

Section 2 Calibrate

Once the freezer has stabilized, the control probe may need to be calibrated. Calibration frequency is dependent on use, ambient conditions and accuracy required. A good laboratory practice would require at least an annual calibration check. On new installations, all parameters should be checked after the stabilization period.

Caution Before making any calibration or adjustments to the unit, it is imperative that all reference instruments be properly calibrated. ▲

Calibrate Control Probe

Plug a type T thermocouple reader into the receptacle located inside the lower door (see Figure 1-5). Compare the control temperature set point to the temperature of the measuring device.

- 1. Press the Mode key until the Calibrate indicator lights.
- 2. Press up/down arrow to match the display to calibrated instrument.
- 3. Press Enter to store calibration.
- 4. Press the Mode key to return to Run mode.

Temperature Stabilization Periods

Startup - Allow 12 hours for the temperature in the cabinet to stabilize before proceeding.

Already Operating - Allow at least 2 hours after the display reaches set point for temperature to stabilize before proceeding.

During calibration, the temperature display will not be available.

If no keys are pressed for approximately five minutes while in calibration mode, the system will reset to Run mode.

Section 3 Alarms

The Model 900 and 5900 Series freezer alarms are displayed on the freezer control panel. When an alarm is active, the indicator next to the alarm description will light and there will be an audible alarm. Press the Silence key to disable the audible alarm for the ringback period. The visual alarm will continue until the freezer returns to a normal condition. The alarms are momentary alarms only. When an alarm condition occurs and then returns to normal, the freezer automatically clears the alarm condition.

Table 3-1. Alarm Indicators

Description	Delay	Ringback	Relay
Power Failure	1 min.	15 min.	Yes
High Temperature Alarm	1 min.	15 min.	Yes
Low Temperature Alarm	1 min.	15 min.	Yes
Probe Failure see 3.2	1 min.	15 min.	No
Door Open	1 min.	15 min.	No
Wrong Power	0 min.	none	Yes
Low Battery*	1 min.	8 hours	No
Low BUS Battery (optional)	1 min.	15 min.	No
Hot Condenser	1 min.	none	No
High Stage Failure	0 min.	15 min.	Yes
Micro Board Failure	0 min.	15 min.	Yes

All alarm delays and ringback times are ±30 seconds.

* The automatic battery test runs immediately on power-up, then every 8 hours thereafter.

High Stage System Failure

The "high stage system failure" condition is created when the high stage compressor and fans run for 30 minutes and are not capable of cooling the interstage heat exchanger to the proper temperature. Under this condition, the high stage compressor and fans will turn off after 30 minutes and an audible and visual alarm will occur. The audible alarm can be silenced and will ring back every 15 minutes.

Probe Failure Alarm

The microprocessor in Model 900 and 5900 Series freezers continually scans all probes including the control probe, heat exchanger probe and condenser probe to ensure that they are operating properly. Should an error be detected, the "Probe Failure" alarm will occur as described above. If an error is detected with the control probe, the high and low stage compressors will run continuously. As a result, the cabinet temperature will decrease until it reaches the lowest temperature that the refrigeration system can maintain. If an error is detected with the heat exchanger probe, the freezer will cycle properly at its temperature set point using a 5 minute step start between the high and low stage compressors. If an error is detected with the condenser probe, there is no impact on the performance of the freezer; however, the hot condenser alarm may also occur. Contact the Technical Services department or your local distributor.

Error	Notes
Er00	Name: Improper model selected.Description: Indicates that DIP SW3 has not selected a proper model or can't be accessed properly.Response: Display shows "Er00" and will not start-up until a proper model is selected. Contact Technical Services.
ErA1	This error condition will prevent peripherals (fans, compressors, etc.) from powering up with the incorrect voltage. Name: Voltage/Frequency failure
ErC1	Description: Indicates the measured RMS line voltage did not agree with the logic level sensed by the micros provided by the installed high voltage PCB; or the measured RMS voltage is not within a tolerable range (180 < Vac < 270 for 230VAC unit / 85 < Vac < 160 for 120VAC unit); or the frequency measured over 10 cycles was not within a tolerable range (55 Hz < Freq < 70 Hz for
Erd1	60 Hz units / 40 Hz < Freq < 55 Hz for 50 Hz units) Response: This condition is checked at power on reset and if it is active, the unit will NOT power up. The unit will indefinitely
ErE1	display "Er_1" in the display and continue to monitor the frequency and voltage. Furthermore, the audible alarm will sound. Other startup error messages may be displayed prior to this message; however, the system will stop the startup sequence for this condition.
ErF1	ErA1 No pulses (zero crossings) detected to determine frequency (50 / 60 Hz) ErC1 Frequency detected is below 50 Hz Erd1 Frequency detected is above 60 Hz (possible noise spikes on supply voltage) ErE1 Unit is 230V and the voltage detected is below the low limit (180VRMS) ErF1 Unit is 230V and the voltage detected is above the high limit (260VRMS)
Erg1	
ErH1	Erg1 Unit is 120V and the voltage detected is below the low limit (85VRMS) ErH1 Unit is 120V and the voltage detected is above the high limit (160VRMS)
Er02	 Name: Control (Cabinet) Sensor Failure Description: This condition indicates that the control sensor has failed to produce a valid reading for ≥12 consecutive reads (~60 seconds). Response: The unit will stage both compressors on (if necessary) and the unit will attempt to head to bottom out. If the sensor recovers, the system will begin to operate normally and respond to the temperature feedback. The remote alarm contacts will become active. 'Er02' will be added to the main display queue and the last valid cabinet temperature value will not be displayed.
Er03	Name: Heat Exchange Sensor Failure Description: This condition indicates that the heat exchange sensor has failed to produce a valid reading for ≥12 consecutive reads (~60 seconds). Response: The display will show "Er03" only when the button sequence to read the heat exchange sensor is depressed.

Error (cont.)	Notes		
Er05	Name: Display Firmware Integrity Failure Description: The display firmware has failed to pass its CRC CCITT checksum integrity test. Response: The display performs this check at startup and the display board will fail to startup without any error indication if it does not pass this at power on.		
Er06	Name: Micro Firmware Integrity Failure Description: The micro firmware has failed to pass its CRC CCITT checksum integrity test. Response: This is checked at power on reset and the "Er06" will be displayed for ~10 seconds at startup if this condition exists.		
Er07	 Name: Micro Fail - CS5521 SPI Failure / UISR Failure Description: This condition indicates a micro board failure due to either the SPI bus is unable to communicate with the ADC device or a UISR event caused the microcontroller to be in an unstable state. Response: The unit will try to recover from this fault three times by a hardware reset of the micro board. In the event that the system couldn't rectify the issue, the following sequence of events will occur: Remote alarm contacts will become active. Buzzer will annunciate audibly and will have a ringback of 15 minutes. "Seven segment" display will show "Er07". The system will have 10 minute staging between the high stage compressor and the low stage compressor activation. The system will go to bottom out temperatures. 		
Er09	Name: Stuck Button Description: This condition indicates that the display board has a stuck button. Response: The Er09 will show on the display periodically.		
Er11	Name: Condenser Probe Sensor Failure Description: This condition indicates that the condenser probe sensor has failed to produce a valid reading for ≥12 consecutive reads (~60 seconds). Response: The display shows "Er11".		
dErr	This is a general display error in which value being displayed can not be represented within characters provided.		
(Four dashes) in display	Name: Lost Communication Description: Communication between the micro board and the display board has been lost. Under this condition, the visual alarm flashes along with dashes in the temperature display (). Contact Technical Services.		

Section 4 Maintenance

^		
[loon	Cabinet	Lytorior
GEAN	Launer	
UIU	U MNIIIUU	

Wipe down the freezer exterior using soap and water and a general use laboratory disinfectant. Rinse thoroughly with clean water and dry with a soft cloth.

Caution Avoid the excessive use of water around the control area due to the risk of electrical shock. Damage to the controls may also result. ▲

Clean Air Filter Th

The air filter should be cleaned four times a year, minimum.

1. Open the front lower door by grasping the bottom left corner.

To remove filter

Figure 4-1. Filter Location

- 2. Locate the grille on the door. See Figure 4-1. Grasp the middle of the grille material and gently pull out to remove.
- 3. Wash the filter material using water and a mild detergent.
- 4. Dry by pressing between two towels.
- 5. Install the filter back into the grille and attach the grille.

Clean Condenser The condenser should be cleaned once per year, minimum.

- 1. Open the front lower door by grasping the bottom left corner.
- 2. Using a vacuum cleaner, exercising care to not damage the condenser fins, clean the condenser.

Depending upon environmental conditions, the condenser may need to be cleaned more frequently.

Clean Water-cooled Condenser	 It is recommended that the water-cooled condenser be cleaned at least once a year, more frequently when environmental conditions are relatively high in particulates. The water-cooled condenser can be cleaned-in-place by using the CIP procedure below. Cleaning solutions can be used, depending on type of deposits or build-up to be removed. Follow manufacturer precautions. Note Do not use liquids corrosive to stainless steel or the brazing material (copper or nickel). Do not use hydrochloric acid or muriatic acid. ▲ 		
CIP (Clean-In-Place) Procedure	1. Disconnect the unit from the water supply.		
	2. Drain the unit.		
	 Drain the unit. Rinse with fresh water and drain the unit again. 		
	4. Fill with fresh water.		
	5. Add cleaning agent (solution and concentration dependent on deposits or build-up).		
	6. Circulate cleaning solution (if feasible).		
	7. Drain the cleaning solution.		
	8. Add and circulate a passivating liquid for corrosion inhibition of plate surfaces.		
	9. Drain this liquid.		
	10. Rinse with fresh water and drain.		
	11. Reconnect the water supply and fill the unit.		
	12. Return to service.		
Prepare Unit for Storage	Defrost the unit as described in 'Defrost Chamber'. This will prepare the unit for storage. Turn off the freezer power switch and disconnect unit from the power source. Turn off the battery power switch (O). Note If the unit has been in service, turn it off and disconnect the power cord connector before proceeding with any maintenance. ▲		

Freezer battery and optional BUS battery

Figure 4-2. Battery Power Switch

Long Term Storage with Water-cooled Condenser



Caution To prevent permanent damage to the refrigeration system of this freezer, proper draining and preparation is required for long term storage. The water-cooled condenser must be dry and sealed. If stored clean and dry, restoring to service only requires a new in-line strainer filter, proper installation and water flow balancing per Section 5, Water Cooled Condenser. ▲

Required: Cleaning solution, supply of dry air or N2

Note Do not use liquids corrosive to stainless steel or the brazing material (copper or nickel). Do not use hydrochloric acid or muriatic acid. ▲

- 1. Turn the freezer off.
- 2. Drain and clean the unit following the procedure in this section, Clean Water-cooled Condenser.
- 3. The condenser must be fully dry prior to storage. Use dry air to purge and dry the condenser circuit as much as possible. Since the freezer is not running, the internal valve will be closed and dry air cannot be circulated through the loop. Allow the condenser to dry, open to the air, until no moisture is seen at the connecting ports.
- 4. After the condenser water loop is dry, purge both the inlet and outlet with nitrogen or dry air and cap the ports.

Note On restarting a freezer that has been stored for some time, it is suggested to run the unit briefly, then check and clean the particulate filter, if fitted, to ensure that no sediment or growth occurred in the condenser during storage. \blacktriangle

Defrost Chamber

- 1. Remove all product and place it in another freezer.
- 2. Turn the unit off and disconnect it from the power source.
- 3. Turn off the battery switch (O). See Figure 4-2.
- 4. Open all of the doors and place towels on the chamber floor.
- 5. Allow the frost to melt and become loose.
- 6. Remove the frost with a soft cloth.
- 7. After defrosting is complete, clean the interior with a non-chloride detergent. Rinse thoroughly with clean water and dry with a soft cloth.
- 8. Plug unit in and turn power switch on.
- 9. Turn the battery power switch to Standby mode (^(U)).
- 10. Allow the freezer to operate empty overnight before reloading the product.

Clean Door Gasket	The door gasket should be cleaned monthly, at minimum. Using a soft sloth, remove any frost build-up from the gasket and door(s). The door gasket may need to be cleaned more frequently if dirt or excessive frost build- up prevents the door from closing properly.	
Vacuum Relief Port	The exterior door gasket provides an excellent seal to protect product, provide an energy efficient thermal barrier to keep cold air in and room temperature air out, and reduce frost build up on the inner doors.	
	Because the door gasket seals so well, a vacuum can be created after a door opening. Warm air enters the cabinet, cools and contracts, creating a vacuum that pulls the door in tightly against the seal.	
	To equalize the pressure inside the cabinet after a door opening requires 1.5- 3.0 cu.ft. of ambient air to be drawn into the cabinet. The amount of air required to equalize the pressure varies depending on the cabinet size, cabinet temperature, duration of door opening, inventory volume and the temperature/humidity of the ambient air. The unit is designed with a "vacuum relief port" that allows the pressure to be equalized.	
	The time required to draw 1.5-3.0 cu.ft. of air into the cabinet depends on	
	two factors: a) size and number of paths available for air to enter the cabinet, and	
	b) pressure difference between internal cabinet and the ambient room.	
	Cabinets with the vacuum relief port operating normally, (i.e. vacuum relief port is not iced over) will require a minimum of 30 seconds up to a maximim of 120 seconds for the cabinet to equalize. This is also a good indication that the exterior door is well sealed.	
	The vacuum relief port requires routine maintenance. It will ice over unless preventive measures are taken. If the vacuum relief port becomes iced over, the freezer will take several hours to equalize pressure.	
	Warning Do not leave the freezer unattended while the door is unlatched. The vacuum could release, resulting in a door opening and product loss. \blacktriangle	
	Observe the inner side of port periodically for frost and ice build-up. Remove any frost with a soft dry cloth. If the tube should become clogged with ice, it must be cleaned. Make sure during cleaning that the vacuum relief tube is completely free of ice to prevent rapid ice formation.	
	Factors that can affect the performance of the vacuum relief port include: high ambient temperature, high humidity conditions and frequent door openings. Maintenance should be performed weekly or as needed.	
	Warning Failure to maintain the vacuum relief port may result in excessive ice build up inside the tube, clogging the port, and inability to open the door. The vacuum relief port may need to be cleaned more often with frequent door openings and high humidity environments.	

Check Battery(s)

- 1. To gain access to the battery, open the lower door by grasping the bottom left corner. The battery is rectangular in shape, located on the front left corner of the compressor compartment and is secured in place by a mounting bracket.
- 2. Directly above the battery(s) is the battery power switch. Turn the battery power switch to the off position (O).
- 3. Cut the tie wrap securing the battery to the mounting bracket. Lift the battery out of the bracket.
- 4. Disconnect the red and black wires from the battery.
- 5. Use a voltmeter set to DC volts. Matching the wire colors, connect the meter to the battery.
- 6. If the voltage reads less than 10.8 volts, replace the battery. If above 10.8, re-install as previously.
- 7. Turn the battery power switch to Standby mode (^(U)).
- 8. Close lower panel door.

Replace Battery(s)

- 1. To gain access to the battery, open the lower door by grasping the bottom left corner. The battery is rectangular in shape, located on the front left corner of the compressor compartment and is secured in place by a mounting bracket.
- 2. Directly above the battery(s) is the battery power switch. Turn the battery power switch to the off position (O).
- 3. Disconnect the battery connections.
- 4. Remove the old battery and install the new battery.
- 5. Reconnect the battery (red to positive and black to negative).
- 6. Turn the battery power switch to Standby mode (^(U)).
- 7. Close lower panel door.

Warning The % of charge can vary depending on the age, usage and condition of the battery. For a consistent and dependable charge, replace the battery every 2 years. Replacement batteries must be rechargeable and are available from Thermo. Refer to the parts list for stock number and description of the replacement batteries. Dispose of the used batteries in a safe manner and in accordance with good environmental practices. ▲

PREVENTIVE MAINTENANCE Freezers

Your equipment has been thoroughly tested and calibrated before shipment. Regular preventive maintenance is important to keep your unit functioning properly. The operator should perform routine cleaning and maintenance on a regular basis. For maximum performance and efficiency, it is recommended that the unit be checked and calibrated periodically by a qualified service technician.

The following is a condensed list of preventive maintenance requirements. See the specified section of the instruction manual for further details.

We have qualified service technicians, using NIST traceable instruments, available in many areas. For more information on Preventive Maintenance or Extended Warranties, please contact the Technical Services Department. Cleaning and calibration adjustment intervals are dependent upon use, environmental conditions and accuracy required.

Tips:

- Fill an upright by starting at the bottom near the probe and add racks to one shelf at a time. Allow freezer to recover to set point between shelves.
- Fill a chest by starting at the left side near the probe. Filling with room temperature racks will result in a long pull-down time.
- Fill unit with frozen product to help overall performance; frozen water jugs, for example.
- Always make certain the vacuum relief port is free of frost and ice, to allow for timely re-entry into the freezer after a door opening.

Action	Monthly	Yearly	Every 2 Years
Verify ambient temperature, <90°F	\checkmark		
* Adjust door handle for firm latching, as needed	v		
Check and clean probe cover, gaskets, hinges and lid(s) of ice and snow. See Figure 1-4 for probe location. See Section 4, Clean Door Gasket.	✓ More frequent cleaning may be required, depending on use and environmental conditions.		
Check air filter. Clean or replace as needed. See Section 4, Clean Air Filter.		✔ 4X	
Check alarm back-up battery. See Section 1, Connect Unit to Electrical Power and Section 4, Replace Battery(s).	<i>v</i>		**Replace
Check condenser fan motor for unusual motor noise or vibration.		~	
* Verify and document calibration, at the minimum, annually. See Section 2, Calibrate.		•	
* Clean condenser compartment and wipe off condenser. See Section 4, Clean Condenser.		~	
Preventive Maintenance *Back-up System - inject test & battery. See Section 5.	v		

To minimize ice build-up inside of freezer:

- Locate the freezer away from drafts or heating/cooling vents
- Keep the number of door openings to a minimum
- Minimize the length of time door is open
- Make sure door latches securely after opening

Section 5 Factory Installed Options

Back Up System (BUS) - P/N 1960445, 1960447

Note Before installation of BUS components, make sure the power to the freezer is disconnected, the battery switch is turned off (O) and the freezer has warmed to ambient temperature. ▲

The built-in BUS (back up system) maintains the freezer chamber temperature below the critical level in the event of a power or equipment failure. If power to the freezer fails, or temperature increases to the back up alarm set point, the BUS injects liquefied gas into the chamber to keep the chamber temperature within the specified range.

The BUS operates on an internal 12-volt, rechargeable battery which is kept charged during normal operation by the integral battery charger.

Install Vent Stack, Solenoid and Injection Assembly

1. Install the injection assembly through the 1/2" pre-punched hole, directly behind the 2" vent stack hole in the center of the chamber ceiling.

Note Cover the open end of injection assembly with tape to keep insulation from entering the nipple. ▲

- 2. Slide 3/8" flatwasher over open end of nipple.
- 3. Insert the covered end of the injection assembly through the exterior hole.
- 4. Remove the tape covering from the end of the nipple and install the 1/8" NPT brass tee on the open end of the nipple. Place Permagum sealant between the brass tee and the interior top.
- 5. Remove the two Phillips head screws securing the metal bracket on the vent stack assembly.
- 6. Install the vent stack through the opening and secure it to the top of the freezer, using screws.

Install Vent Stack, Solenoid and Injection (continued)

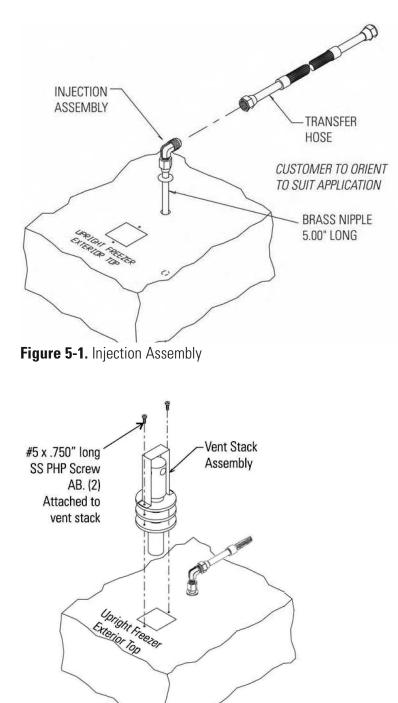


Figure 5-2. Vent Stack Assembly


Install Vent Stack, Solenoid and Injection (continued)

Install Temperature Probe

- 7. Inside the chamber, seal around the end of the vent stack with Permagum.
- 8. Install the transfer hose connecting one end to the injection assembly, the other to the solenoid valve. Connect the solenoid valve to the supply source. The solenoid mounting bracket may be discarded.

Note When selecting a CO₂ supply cylinder, it must be equipped with a siphon tube. \blacktriangle

- 1. Locate the 0.500" pre-punched hole in the upper left hand back corner of the chamber ceiling. Remove the tie wrap securing the coiled probe/solenoid harness. Uncoil the probe lead and run the probe tip (approximately 12") down through 0.500" porthole (Figure 5-4).
- 2. As shown in Figure 5-3, thread the small tie wrap through the openings in the front of the bracket. Secure the probe on the back of the bracket with the tie wrap.
- Tap #8-32 the two pre-punched holes located on the interior left wall of the freezer. Mount the bracket. Figure 5-4 shows the Back-Up probe mounted on the interior left side wall of the freezer.

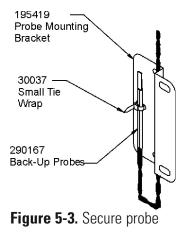

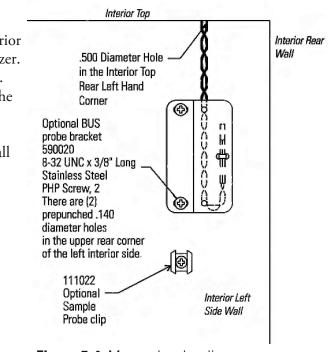


Figure 5-4. Mounted probe clip

2

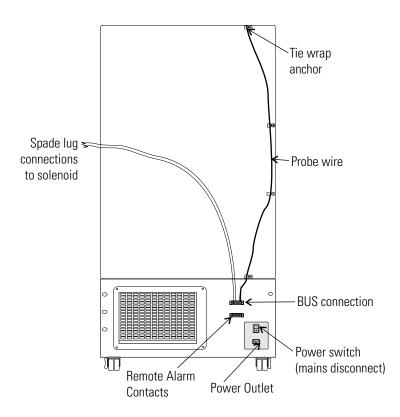


Figure 5-5. Probe wire and solenoid connections

1. Remove the four screws on the freezer back panel and use them to mount the tie wrap anchors as shown in Figure 5-5. Secure the probe wire with tie wraps.

- 2. Plug the solenoid/probe connector into the BUS connection and secure with a screw on the right and left side. The connector is keyed.
- 3. Loosen the terminal screws on the solenoid. Slide the spade lug connectors under the screws and tighten to secure.
- 4. Connect power to the freezer. Turn the freezer On, with battery switch Off (O).

a. Arming on Power On

The arming function occurs only one time at initial power-on. It may occur in 3 ways; BUS timeout arming, both sensors fail and BUS temperature sensor arming.

At power-on, the BUS will perform a Hold-off period check. The BUS will not inject the refrigerant for a Hold-off period, calculated as follows:

Connect Probe/Solenoid Harness

BUS Control Panel

Hold-off period:

- 12 hour is between 0 and +70C
- 12 hour is between 0 and (BUS Setpoint Hysteresis)
- 0 hours is below (BUS Setpoint Hysteresis)
- 0 hours is temperature Out Of Range or both probes fail.
- b. The Low Battery indicator may also illuminate.
- 5. Turn the battery switch to Standby mode (^(U)) to charge both batteries.

The following section describes the configuration and operation of the BUS.



Warning When activated, this unit injects liquid nitrogen or carbon dioxide. Liquid nitrogen can cause serious freezing (frostbite) if it comes in contact with unprotected skin or eyes. Nitrogen suppresses oxygen levels and may cause suffocation if area is not well ventilated. Refer to Appendix A for the proper handling of liquid LN₂.



Caution Make sure the pressure relief valve on any LN₂ tank is adjusted to 30 PSI maximum blow-off. ▲



Warning Carbon dioxide gas suppresses oxygen levels and may cause suffocation if area is not well ventilated. Refer to "Handling Liquid CO₂ in Appendix B of this manual. ▲

Backup System



Figure 5-6. BUS Control Panel

Power - indicates the unit has AC power.

- 1. Normal Operation
- Low Battery: Illuminates only after a battery test which occurs every 8 hours. The battery test is a loaded test and during the test, the battery voltage is monitored. When the battery voltage is less than a certain threshold, this light will illuminate and stay illuminated until the next battery test.
- Solenoid Engaged: Illuminates any time the controller is activating the solenoid. Illumination will only occur when the controller senses the solenoid is active and calling for injection.

BUS Control Panel (cont.)

2. Fault Code

This same user interface will also permit fault codes to be easily viewed when the user presses and releases the Press to Test button in less than 2 seconds. When this occurs, both the Low Battery LED and Solenoid Engaged LED will blink twice ON/OFF quickly indicating that fault code display mode is active. For the next ~10 seconds, the fault codes will be displayed as shown in the table below. The display shows only the lower number fault code until that fault is resolved.

Low Battery LED	Solenoid Engaged LED	Fault Code Number	Fault Code Definition
OFF	OFF	0	Normal Operation
OFF	Blinks 1 time	1	Primary Probe Failure
OFF	Blinks 2 times	2	Secondary Probe Failure
OFF	Blinks 3 times	3	Solenoid Driver Failure
OFF	Blinks 4 times	4	Low Tank

Table 5-1. BUS Fault Codes

Each blink will take one second, with half of the second interval being ON (illuminated) and the other half of the second interval being OFF. After the sequence concludes, both the Low Battery LED and the Solenoid Engaged LED will flash four times ON/OFF quickly to indicate a return to normal display.

3. BUS Injection Test

The user must press and hold the Press-to-Test button for ≥ 2 seconds before injection occurs. After the Press-to-Test button is pressed for ≥ 30 seconds, the solenoid will turn off. This will avoid a button or button signal becoming stuck in an active condition.

During a Press-to-Test injection, if either RTD #1 or RTD #2 reads a temperature less than 20^oC below BUS set point and it is determined to be a valid input, then injection will terminate immediately.

4. BUS Set point "Validation" Feedback

Whenever the BUS set point or working fluid changes, the BUS acknowledges the set point acceptance by causing the Low Battery LED and Solenoid Engaged LED to blink three times ON/OFF quickly. This does not imply that the BUS error-checked the values, other than that the value received was within the range of the BUS working correctly with prior programming.

Note Solenoid will not engage if door is open. ▲

BUS Control Panel (cont.)

Caution The back-up system is designed to inject refrigerant (CO₂ or LN₂). In the unlikely event of both probes failing, the back-up system will operate in PWM or Non-PWM mode based on jumper setting (Figure 5-7). \blacktriangle

<u>!</u>

For Hysteresis & Dual probe fault mode settings (Figure 5-7), locate Jumper J17 near to solenoid/ sensor connector (J1), refer to table below for operation details:

Jumper PIN	Status	Hysteresis	Dual Probe Fault Mode
J17 PIN 1 & PIN 2	Short Together/ Jumper Cap Installed	5°C	Non-PWM Mode
J17 PIN 1 & PIN 2	Open/No-Jumper Installed	9°C	PWM Mode

 Table 5-2.
 Probe Fault Modes

<u>PWM mode</u>: If both probes fail, the solenoid will activate peridically, allowing the flow of refrigerant (CO2 or LN2).

<u>Non-PWM mode:</u> If both probes fail, the solenoid remains OFF and no refrigerant will flow.

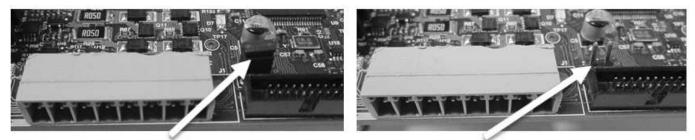


Figure 5-7. Jumper Cap for Hysteresis & Dual Probe Fault Mode Selection

Configure Optional BUS (Back-Up System)

The optional BUS can be configured for LN_2 or CO_2 supply. To select the supply type:

- 1. Press the Mode key until the Backup indicator lights.
- 2. Press the up or down arrow key. The display will show "cO2" or "Ln2" for selection.
- 3. Press Enter to save the setting.
- 4. Press the Mode key until the Run indicator lights for Run mode

If no control keys are pressed, the freezer will automatically return to to RUN mode after 5 minutes.

Set Optional BUS Set Point

The optional back up system is designed to inject CO_2 or LN_2 into the freezer compartment if the temperature rises above back up system set point. To set the BUS set point:

- 1. Press the Mode key until the Set Temperature and Backup indicators light.
- 2. Press the up or down arrow key until the desired BUS set point is displayed.
- 3. Press Enter to save the setting.
- 4. Press the Mode key until the Run indicator lights for Run mode

If no control keys are pressed, the freezer will automatically return to to RUN mode after 5 minutes.



Warning Changing the operating temperature set point can affect the BUS set point. The BUS set point will self-adjust to maintain a temperature of at least 10°C above the operating temperature set point. ▲

Note The BUS set point can not be set any colder than the high temperature alarm set point. (See Section 1, Set High Temperature Alarm). If the back-up system is installed with CO_2 , then -65°C is the coldest BUS set point that can be used (if the cabinet set point is -75°C or colder).

After the freezer has stabilized and both batteries are fully charged, the BUS can be tested to verify proper operation.

Test BUS Operation

- 1. Disconnect the AC power to the freezer by turning power switch off.
- 2. As the freezer warms up, verify the BUS injects at the desired temperature. Displayed temperature may vary by a few degrees from inject temperature due to the differences in probe locations.
- 3. <u>Preventive Maintenance</u> Monthly maintenance action to check CO₂ or LN₂ back-up system operation, and battery voltage.
 - Check for fault codes on the BUS probe per Table 5-1. Also verify solenoid operation by performing a Press-To-Test.

Clean Vent Stack

Routinely check the vent stack for frost or ice build-up. The type of frost that forms in the vent stack is generally very soft and may be easily removed with a bristle brush or soft cloth. If ice build-up has occurred, a complete defrost may occasionally be required. See Section 4 for freezer defrost instructions.

Disconnect Fitting Asm. and Transfer Hose

- To disconnect the freezer back-up from the gas supply:
- 1. Close the supply valve.
- 2. Depress the test button on the BUS control box to remove the gas from the line.
- 3. <u>Slowly</u> disconnect the fitting assembly from the supply (in the event that any gas remains in the line).

Chart Recorder

The following section describes the set up and operation of the optional chart recorder.

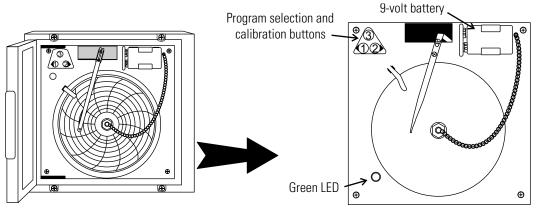


Figure 5-8. Recorder Details

Installing the Chart Paper

- 1. Open the plastic door of the recorder and press button #3 until the pen begins to move outward.
- 2. Unscrew the knob at the center of the chart and remove the paper.
- 3. Install the new chart paper, position the paper to the correct time line and replace the knob.
- 4. Remove the cap from the felt pen and press button #3.

Change Recorder Temperature Range

The chart recorder contains eight factory programmed temperature ranges. To change the recorder range:

- 1. Press and hold button #3 for one second, then let the pen move off the chart paper.
- 2. Press and hold for five seconds either button #1 or button #2.

Change Recorder		3. Release button and the green LED begins to		From	To
Temperature Range	flash. Count the number of flashes to		1	-40	30°C
(cont.)	determine the present p	determine the present program setting.		0	60°C
	4. To change the program		3	-100	38°C
	or right arrows to incre	ase or decrease the	4	-5	50°C
		count.		0	100°C
	5. When the desired progr		6	-100	200°C
	flashing, press button # arm back onto the char	e 1	7	-115	50°C
	begin in the new progra	e	8 Table 5 3	-10	70°C
			Table 5-3 Range	D. NECUN	uer remp
Recorder Calibration	Note The recorder must be following calibration proceed		rs before p	perform	ing the
	1. Place an accurate thermometer in the chamber next to the recorder probe. Temperature probes for the recorder are located in the left front corner of the freezer chamber (Figure 1-4).				
	2. After about three minutes, compare the thermometer reading with the chart recorder reading.				
	3. If an adjustment is necessary, press the #1 button to move the pen to the left or the #2 to move the pen to the right. The button must be held about five seconds before the pen begins to move. Release the button when the pen position matches the thermometer.				
	Note The felt-tip pen on the recorder requires periodic replacement. The ink will appear to fade before replacement becomes necessary. Additional pen tips may be purchased from Thermo. ▲				
Water-cooled Condenser	The water-cooled condenser is a factory installed option (195964 [13 c ft], 195965 [17, 23, 28 cu ft], 195967 [12, 17, 20 cu ft chest]) and requires a qualified technician at freezer installation. Refer to Table 5-4 the specifications for this option.			d	
	Table 5-4. Water-Cooled Content	idenser Specifications			
	Water Pressure	Not to exceed 90 ps	si (620.5 kpa)	
	Water Temperature Range	Water Temperature RangeNot to exceed 29.4°C (85°F)			
	Inlet Connection	Inlet Connection 1/2" NPT			
	Outlet Connection	1/2" NPT			

Flow Rate Required (minimum)

1.0 gallons (3.8 liters) per minute

Water-cooled Condenser (cont.)

Water Quality: Free of particulates that could cause a blockage, or impair function of the regulating valve or heat exchanger. A stainless steel inline strainer is factory installed in the inlet pipe to minimize particulates in the water supply.

Water Drainage: If water return line labeled WATER OUTLET is not connected to a closed loop chiller system, connect the return line water to a reservoir tank or drain. This will hold/dispose of the water from the freezer to avoid flooding.

To verify operating requirements, you will need:

- a <u>flowmeter</u> capable of measuring up to at least 5.3 gallons (20 liters) per minute to measure flow rates at the water return line labeled WATER OUTLET.
- a <u>pressure gauge</u> ranging from 0 to 300 psig to measure pressure at the water inlet labeled WATER INLET.
- a <u>thermocouple or other temperature monitor</u> to measure the water temperature near the water inlet fitting labeled WATER INLET.

Clearance: Before installing the water connections, position the freezer with an extra 2 inches (5cm) clearance in back for a hose. This is in addition to the 6 inches (15cm) clearance in back previously recommended in this manual.



Caution When using multiple freezers in a water loop, proper water flow and inlet temperature must be maintained for EACH freezer in the loop. It is recommended that a valve be installed in the supply line of each freezer to facilitate a balanced flow rate.


If the number of freezers in a loop is more than 5, it is recommended that one freezer be started at a time and tested to verify water flow, before starting all and allowing to run. Water flow of at least 1 gpm, at or below the maximum allowed inlet temperature, is required. Values higher than minimum are acceptable. Do not allow the freezer to continue to run with no observed water flow within several seconds of the compressor starting. ▲


Recommended Bypass Water Loop Assembly: The bypass loop contains fittings, a valve, and a flowmeter to simulate the presence of a freezer and allows for balancing of the entire water supply system to ensure that 1 gpm is available for each freezer in the loop before the freezers are started. Values higher than minimum are acceptable. **Do not adjust the valve inside the deck of the freezer. It is pre-set at the factory.**

Water-cooled Condenser (cont.)

Option A: Install the bypass instead of a freezer (Figure 5-9).

<u>Option B</u>: Install the bypass permanently (Figure 5-9) so that all freezers can be left connected but water flow can be simulated using the bypass without turning on all freezers.

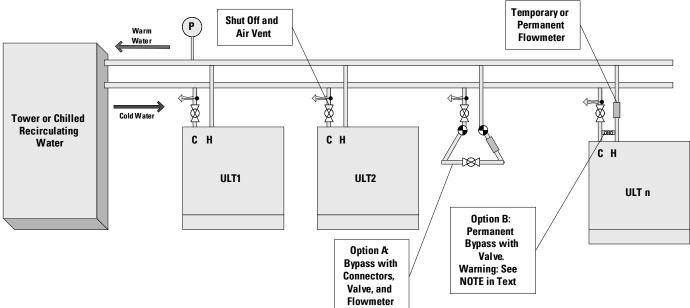

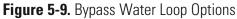



Water Connections

Caution If Option B is used, the bypass valve MUST be closed prior to freezer operation to ensure all water passes through the freezer rather than the bypass. Failure to do this will result in damage to the freezer.

To install the water connections:

- 1. It is recommended to install automatic air purge vents near the water inlet (or in each branch circuit in larger systems) to ensure air is purged from the lines as the system operates.
- 2. Remove the caps from the WATER INLET/WATER OUTLET fittings on lower back of the freezer.

Caution Thread sealant or Teflon^{\circ} tape is highly recommended for use on threaded fittings. \blacktriangle

- 3. Make sure the water line is clear of blockage or debris.
- 4. Connect the incoming water supply line to the WATER INLET fitting.
- 5. Connect the water return/drain line to the WATER OUTLET fitting.
- 6. Adjust to specified operating standards (pressures, temperature).



Water Connections (cont.)

- 7. Turn on water flow to the unit. Check for leakage. Close bypass, if installed.
- 8. Refer to the Installation and Start-Up section of this manual **before** starting the freezer.
- 9. Adjust the water flow rate. Note that the water will not flow through the outlet until the first stage compressor is running. Do not confuse the fan running with the compressor start-up.



Figure 5-10. Inlet/Outlet Fittings



Caution Only authorized service

personnel should turn the freezer on without water connected and flowing.

Caution See Bypass Water Loop Assembly information, if required. \blacktriangle

Five Inner Door Option

The five inner door option (P/N 189405 [13 cu ft], 189406 [17 cu ft], 189407 [23 cu ft], 195652 [28 cu ft]) is factory installed. The freezer is converted to accommodate four adjustable specimen shelves with the fifth "shelf" as the bottom of the freezer chamber.

Section 6 Specifications

Single Door Units

Catalog No.	902	903	909	
Temperature Range	-50°C (-58°F) to -86°C (-123°F) in an 18C to 32C * (64.4F to 89.6F) ambient			
Exterior Dimensions	33.3"W x 77.9" H x 32.9" F-B (84.6cm x 197.9cm x 84.6cm)			
Interior Dimensions	23.0"W x 51.5"H x 19.3" F-B (58.4cm x 130.8cm x 49.0cm)			
Capacity	13.0 cu. ft. (368.1 liters)			
Refrigeration	Two 1 HP (2545 BTUH each)			
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door			
Electrical - nominal voltage ±10%	230VAC, 50 Hz, 12.0 FLA Operating Range: 208VAC-240VAC	120VAC, 60 Hz, 16.0 FLA Operating Range: 108VAC-130VAC	208-230VAC, 60 Hz, 12.0 FLA Operating Range: 208VAC-240VAC	
Breaker Requirements	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	20 Amp, Dedicated Circuit, 20 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	
Shipping Weight	696 lbs. (316 kg)			

Catalog No.	904	905	910	
Temperature Range	-50°C(-58°F) to -86°C(-123°F) in an 18C to 32C * (64.4F to 89.6F) ambient			
Exterior Dimensions	33.3"W x 77.9" H x 38.9" F-B (84.6cm x 197.9cm x 94.0cm)			
Interior Dimensions	23.0"W x 51.5"H x 25.3" F-B (58.4cm x 130.8cm x 64.3cm)			
Capacity	17.3 cu. ft. (489.9 liters)			
Refrigeration	Two 1 HP (2545 BTUH each)			
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door			
Electrical- nominal voltage ±10%	120VAC, 60 Hz, 16.0 FLA Operating Range: 108VAC-130VAC	230VAC, 50 Hz, 12.0 FLA Operating Range: 208VAC-240VAC	208-230VAC, 60 Hz, 12.0 FLA Operating Range: 208VAC-240VAC	
Breaker Requirements	20 Amp, Dedicated Circuit, 20 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	
Shipping Weight	763 lbs. (346 kg)			

Section 6 Specifications

Cotolog No	006		OEC		057		
Catalog No.	906		956		957		
Temperature Range	-50°C (-58°F) to -86°C (-123°F) in an 18C to 32C * (64.4F to 89.6F) ambient						
Exterior Dimensions	40.8"W x 77.9" H x 38.9" F-B (103.6cm x 197.9cm x 94.0cm)						
Interior Dimensions	30.6"W x 51.5"H x 25.3" F-B (77.7cm x 130.8cm x 64.3cm)						
Capacity	23.0 cu. ft. (651.3 liters)						
Refrigeration	Two 1 HP (2545 BTUH each)						
Insulation	Non	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door					
Electrical- nominal voltage ±10%	230VAC, 50 Hz, 12.0 FLA Operating Range: 208VAC-240VAC		120VAC, 60 Hz, 16.0 FLA Operating Range: 108VAC-130VAC		208-230VAC, 60 Hz, 12.0 FLA Operating Range: 208VAC-240VAC		
Breaker Requirements	 15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker 		20 Amp, Dedicated Circuit, 20 Amp Time Delay Breaker		15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker		
Shipping Weight	900	lbs. (408 kg)					
Catalog No.		907		908]	
Temperature Range	-50°C (-58°F) to -86°C (-123°F) in an f			an 18C to 32C * (64.4F to 89.6F) ambient			
Exterior Dimensions		46.8"W x 77.9" H x 38.9" F-B (118.9cm x 197.9cm x 94.0cm)					
Interior Dimensions	36.6"W x 51.5"H x 27.0" F-B (93.0cm x 130.8cm x 68.6cm)					-	
Capacity	28.0 cu. ft. (792,8 liters)					-	
Refrigeration	Two 1 HP (2545 BTUH each)					-	
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door						
Electrical- nominal voltage ±10%		230VAC, 50 Hz, 12.0 FLA Operating Range: 208VAC-240VA		208-230VAC, 60 Hz, 12.0 FLA AC Operating Range: 208VAC-240VAC		-	
Breaker Requirements		15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker		15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker			
Shipping Weight		1070 lbs. (485 kg)					
Catalog No.	990		995		989		
Temperature Range	-50°	C(-58°F) to -86°C(-123°F) in an 18	BC to	32C * (64.4F to 89.6F) ambient			
Exterior Dimensions	40.8"W x 77.9" H x 38.9" F-B (103.6cm x 197.9cm x 94.0cm)						
Interior Dimensions	30.6	"W x 51.5"H x 25.3" F-B (77.7cm	x 130).8cm x 64.3cm)			
Capacity	23.0 cu. ft. (651.3 liters)						
Refrigeration	Two 1 HP (2545 BTUH each)						
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door						
Electrical - nominal voltage ±10%		/AC, 60 Hz, 16.0 FLA rating Range: 108VAC-130VAC	230VAC, 50 Hz, 12.0 FLA208-230VAC, 60 Hz, 12.0 FLAOperating Range: 208VAC-230VACOperating Range: 208VAC-230VAC				
Breaker Requirements	20 A	mp, Dedicated Circuit, mp Time Delay Breaker			15 Amp,	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	
Shipping Weight	900	lbs. (408 kg)	1		1		

Double Door Units

Catalog No.	991	992	986				
Temperature Range	-50°C (-58°F) to -86°C (-123°F) in an 18C to 32C * (64.4F to 89.6F) ambient						
Exterior Dimensions	33.3"W x 77.9" H x 32.9" F-B (84.6cm x 197.9cm x 83.6cm)						
Interior Dimensions	23.0"W x 51.5"H x 19.3" F-B (58.4cm x 130.8cm x 49.0cm)						
Capacity	13.0 cu. ft. (368.1 liters)						
Refrigeration	Two 1 HP (2545 BTUH each)						
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door						
Electrical- nominal voltage ±10%	230VAC, 50 Hz, 12.0 FLA Operating Range: 208VAC-240VAC	120VAC, 60 Hz, 16.0 FLA Operating Range: 108VAC-130VAC	208-230VAC, 60 Hz, 12.0 FLA Operating Range: 208VAC-240VAC				
Breaker Requirements	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	20 Amp, Dedicated Circuit, 20 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker				
Shipping Weight	696 lbs. (316 kg)						
Catalog No.	993	994	988				
Temperature Range	-50°C (-58°F) to -86°C (-123°F) in an 18C to 32C * (64.4F to 89.6F) ambient						
Exterior Dimensions	33.3"W x 77.9" H x 38.9" F-B (84.6cm x 197.9cm x 94.0cm)						
Interior Dimensions	23.0"W x 51.5"H x 25.3" F-B (84.6cm x 130.8cm x 64.3cm)						
Capacity	17.3 cu. ft. (489.9 liters)						
Refrigeration	Two 1 HP (2545 BTUH each)						
Insulation	Non CFC, foamed-in-place urethane: 5.0" (12.7cm) cabinet; 4.5" (11.4 cm) door						

230VAC, 50 Hz, 12.0 FLA

	voltage ±10%	Operating Range: 108VAC-130VAC	Operating Range: 208VAC-240VAC	Operating Range: 208VAC-240VAC	
-	Breaker Requirements	20 Amp, Dedicated Circuit, 20 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	15 Amp, Dedicated Circuit, 15 Amp Time Delay Breaker	
	Shipping Weight	763 lbs. (346 kg)			
-					

120VAC, 60 Hz, 16.0 FLA

Electrical- nominal

208-230VAC, 60 Hz, 12.0 FLA

Section 6 Specifications

Certifications

Declaration of Conformity is available upon request

Safety Specifications

Indoor Use Only Altitude - up to 2,000 meters Temperature - 5°C to 43°C Humidity - Maximum RH 80% for temperatures up to 31°C, decreasing linearly to 50% RH at 40°C Mains Supply Fluctuations - Mains supply voltage fluctuations not to exceed ±10% of the nominal voltage Installation Category II ¹ Pollution Degree 2² Class of Equipment I

Intended Use

This product is intended for use as a General Purpose Laboratory Freezer for storing samples or inventory between -50 and -86C.

Unintended Use

This unit is not intended for use in an explosive environment, nor to be used for the storage of flammable inventory. This unit is not intended for use in a Class II medical application as defined by Title 21 of the Federal Code of Regulations.

Fluorinated Greenhouse Gases

Compliant with Regulation(EU) No 517/2014 of the European Parliament and of the council on fluorinated greenhouse gases.

This product contains fluorinated greenhouse gases in a hermetically sealed system. This product contains foam blown with fluorinated greenhouse gas, R-245fa. If a leak in the sealed system is detected, the operator shall repair without undue delay.

For ULT models rated 230V/50Hz, intended for import into the EU, the following sizes are designed for use with the following amounts of fluorinated greenhouse gases:

		1st Stage		2nd Stage		
Size (ft ³)	Refrigerant	Amount (kg)	GWP	Refrigerant	Amount (kg)	GWP
13	R-404a	0.723	3 922	R-508b	0.306	13 396
17	R-404a	0.751	3 922	R-508b	0.388	13 396
23	R-404a	0.751	3 922	R-508b	0.425	13 396
28	R-404a	0.652	3 922	R-508b	0.397	13 396

¹ Installation category (overvoltage category) defines the level of transient overvoltage which the instrument is designed to withstand safely. It depends on the nature of the electricity supply and its overvoltage protection means. For example, in CAT II which is the category used for instruments in installations supplied from a supply comparable to public mains such as hospital and research laboratories and most industrial laboratories, the expected transient overvoltage is 2500V for a 230V supply and 1500V for a 120V supply.

² Pollution degree describes the amount of conductive pollution present in the operating environment. Pollution degree 2 assumes that normally only non-conductive pollution such as dust occurs with the exception of occasional conductivity caused by condensation.

thermoscientific.com

© 2014 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

Thermo Fisher Scientific (Asheville) LLC 401 Millcreek Road Marietta, Ohio 45750 United States

